

PneuWalker™

MOBILITY TRAINER



The PneuWalker Mobility Trainer is the latest innovation in mobility technology from Pneumex. It offers freedom from the fear of falling for both the patient and the therapist. On the PneuWalker the patient can experience natural movement while working vertically on proper gait patterns.



Patent Pending

EXTRAORDINARY VERSATILITY

- User-friendly, analog control allows for precise unweighting.
- Height adjustable from 3' 6" to 6' 6".
- Width adjustable to go through doorways or over treadmills.
- Unweights unilaterally to address hemiplegic conditions.
- Onboard pneumatic compressor with rechargeable battery that runs up to 12 hours on a single charge.
- Adjustable handlebars swing out of the way for easy access and patient transfer.
- Allows a solo clinician to lift and transfer patients.
- Patented X-Heart Vest with adjustability for posture and gait.

A solo clinician can safely lift, transfer and train their patient

360 Degree Swivel attachment allows for 360 degree rotation



Key Benefits

- Abilities patients
- Bariatric patients
- Orthopedic patients
- Strokes
- Back problems
- Gait training
- Cardiovascular
- Respiratory ailments
- Neurological patients
- Safety



Analog Control Panel

Easy to read, easy to use analog control allows precise unweighting from 0 to 300 lbs.



PneuWalker Maximum lift 300 pounds
Pneumex recommends a maximum patient weight 400 pounds

"The PneuWalker, with the optional PneuGait Straps and Gait Assist Bands, allow me to treat my hemiplegia patient instead of spending all of my efforts stabilizing him." Kim, PT

"This is the best combination vest and unweighting system we have found for the long term care." Dusty, PTA

"The PneuWalker is the first system I have been on in 26 years that has not caused dysflexia Brian (Partial Quad)

Unweighting Studies on the Web

www.pneumex.com

Genuine Pneumex Equipment

www.pneumex.com

pneumex@pneumex.com

800.447.5792

Specifications

ACCESSORIES

X-Heart PneuVest

A totally new concept in unweighting. Light weight, comfortable, color coded, and easy to use. Securely holds and lifts patients with the primary lift point at the posterior rib cage at T12.



One-PneuGait Accessory Pack

- 1-Single Ankle Cuff
- 1-Double Ankle Cuff (for DropFoot)
- 1-PneuGait Foot Strap with Buckle
- 1-Gait Assist Band



Safety Drop Stop

Limits the patient drop to 2" if they stumble and start to fall.

Rechargeable Power Supply

Overnight charge provides sufficient power to operate all day, up to 12 hours.

Lifting Bar


Uniquely designed to unweight more on one side than the other which can be beneficial when dealing with patients who are primarily affected on one side.

PneuWalker 360 Degree Swivel Attachment

Attaches to PneuWalker to allow for 360 degree rotation



SPECIFICATIONS::

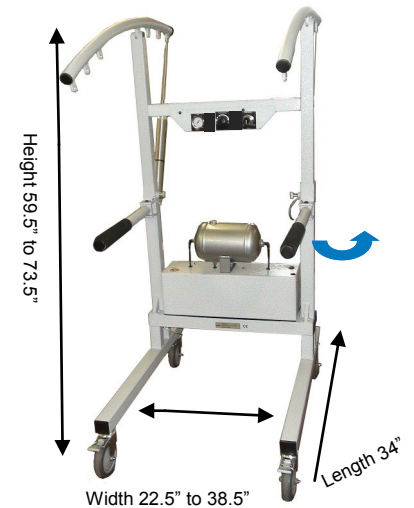

ADJUSTABLE FEATURES:



Adjustable Width

Easily expand PneuWalker inside width to fit over treadmills and balance equipment. An easy adjustment reduces width to go through virtually any doorway.



Adjustable Handlebars

Swing out bars allow for easy transfers from wheelchair, bed or commode.



Height		Width Outside		Width Inside	
Min	Max	Min	Max	Min	Max
59.5 in	73.5 in	26.5 in	42.5 in	22.5 in	38.5 in
151.1 cm	186.7 cm	67.3 cm	108 cm	55.9 cm	97.8 cm

Lifting	Range	Handle Bars	Clearance
Max		2 in Adj	Floor up
300 lb	20-29 in	30-40 in	12 in
136.36 kgs	50.8-73.7 cm	76.2-101.6 cm	30.48 cm

Distributed by:

Pneumex
2605 North Boyer Ave
Sandpoint, Idaho 83864
800.447.5792



Genuine Pneumex Equipment

www.pneumex.com

pneumex@pneumex.com

800.447.5792