

# ULTRA PRO™ Ultrasound Table

Comfort-focused design for  
quality imaging



## Features

- Ergonomic design for enhanced patient and sonographer comfort
- 500 lb patient capacity, suitable for bariatric patients
- Drop-down leg section for chair-like comfort
- Height range down to 23" accommodates wheelchair transfers and those with ambulatory difficulties
- Fowler positioning, 0 to 80 degrees, designed to accommodate natural body extension and sitting position without slide or shift
- Motorized height, Fowler positioning, Trendelenburg motions and auto-level adjustment
- Naugahyde, antimicrobial mattress cover with Advanced Beauty Gard provides protection against bacteria
- Central floor-locking system maintains table stability
- Durable hand controller activates table positions
- Allows seated or standing work positions

## Optional

- Foot Controller is ideal for when hands are busy with other functions
- Side rails mount flush to table for unencumbered patient access
- Stirrups are designed to cup at the heel for added comfort
- Articulating Scanning Arm Board, adjustable from 0 to 130 degrees
- Headrest for easier thyroid and carotid access
- IV Pole easily attaches for quick convenience
- Paper Dispenser accommodates a roll of hygienic table paper. Includes a cutter strap.

## Description

The Ultra Pro Ultrasound Table is the professional's choice when performing general ultrasound procedures, including OB/GYN.

The table's ergonomics are designed to reduce sonographer fatigue and minimize patient repositioning, making it easy to achieve quality images in less time.

Fowler positioning, which provides the utmost in patient comfort, an enhanced central floor-locking system offering easier access from around the table, and an antimicrobial vinyl mattress ensuring safety against infection are all standard on the Ultra Pro table.

Customization is made possible with optional items including flush-mounted side rails allowing complete patient access and safety, retractable stirrups with contoured heel cups and an extra-long extension for tall patients, an articulating scanning arm board and a headrest that drops down 15 degrees to allow for carotid studies.

## Biodex | Ultra Pro Ultrasound Table



All Ultrasound Tables feature an adjustable height range of 23" to 39" – ideal for wheelchair transfers and patients with ambulatory difficulties.

### Specifications

#### Dimensions:

- **Overall:** 70" l x 30" w (177.8 x 76.2 cm); 35" w (88.9 cm) with optional side rails
- **Weight:** 335 lb (152 kg)
- **Drop-Down Leg Section:** 12.6" l (32 cm)

#### Motions:

- **Height Adjustable:** 23" to 39" (58.4 to 99 cm)
- **Trendelenburg:** 0° to +15°, 0° to -15°
- **Fowler Positioning:** 0° to 80° infinitely adjustable

#### Controls:

- **Hand Control:** Activates height, Trendelenburg motions, Fowler positioning and auto level motions
- **Foot Control (optional):** Activates height, Trendelenburg motions and Fowler positioning
- **Tabletop:** Three primary sections: torso section with Fowler positioning, center section remains fixed and leg section drops down to 40° and 80° for stirrup access
- **Wheels:** 5" (12.7 cm) Individual locking swivel casters; central floor-locking system
- **Patient Restraints:** Two body straps
- **Mattress:** Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; Naugahyde, antimicrobial mattress cover with Advanced Beauty Gard provides protection against bacteria
- **Upholstery Color:** Graphite
- **Finish:** Powder coat
- **Patient Capacity:** 500 lb (227 kg); weight tested to four times the patient load rating
- **Power:** 115 VAC or 230 VAC



Optional flush-mounted side rails provide added safety without obstructing patient access.

- **Warranty:** Two years parts and labor
- **Certifications:** ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2 No. 601.1.M90 and EN 60601-1 standards
- All ultrasound tables are compatible with Hoyer Lifts
- **Primary Function:** General Ultrasound
- CDRH registered and FDA Device Listed

### Ordering Information


- 058-720 Table, Ultrasound, Ultra Pro, 115 VAC  
Includes hand controller.
- 058-725 Table, Ultrasound, Ultra Pro, 230 VAC  
Includes hand controller.

### Related:

- 058-704 Chair, Ergonomic, Sonography

### Optional:

- 058-741 Foot Controller
- 058-633 Side Rails, Folding
- 058-652 Stirrups, Retractable
- 058-736 Arm Board, Articulating Scanning
- 058-737 IV Pole
- 058-611 Dispenser, Paper
- 058-738 Headrest

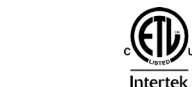

[www.biodex.com](http://www.biodex.com)